


## Setting Up the SmartMonitor 2

### Step 1: Set the monitor on a clean flat surface and gather you supplies and equipment.

- Gather the following supplies/equipment: Smart Monitor, lead wires, electrodes, belt
- Be sure the speaker is not blocked.
- Makes sure the monitor is close enough to connect to the baby comfortably.


### Step 2: Connect the patient cable to the monitor.

- Insert the round end of the patient cable into the bottom round connector found on the front of the monitor.
- Line up the notch on the connector and push until you feel the connector snap into place.
- To remove the patient cable, grasp the patient cable at the base of the patient input connector and gently pull back. Do not grasp and pull the patient cable to remove from the monitor.

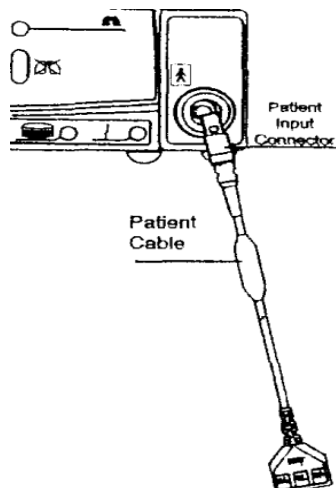


### Step 3: Connect the lead wires to the patient cable.

- The larger end of the patient cable has 3 openings, marked **LA** (black), **RL** (green), and **RA** (white).
- Take the **white** lead wire and insert it into the opening marked **RA**.
- Take the **black** lead wire and insert it into the opening marked **LA**.
- Firmly push each lead wire in until the wire snaps into place.

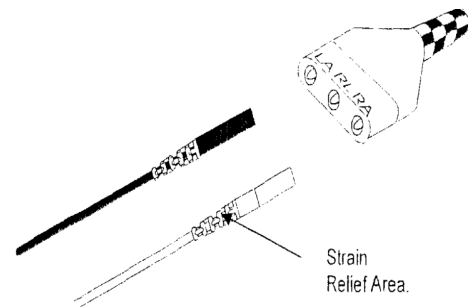

### Step 4: Connect the lead wires to the electrodes.

- Insert the black **LA** lead wire into one electrode.
- Insert the white **RA** lead wire into the other electrode.
- Make sure the metal tips of the lead wires are fully inserted into the electrodes.

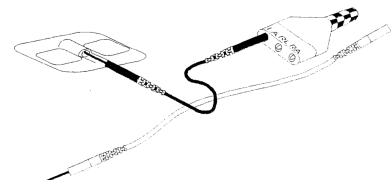


If you have any questions or concerns, please call Cook Children's Home Health at (800) 747-8242. A nurse is available 24 hours a day.

## Setting Up the SmartMonitor 2


### Step 5: Attach the electrodes to the infant belt.

- Place the electrode belt on a flat surface.
- Lay your baby on the belt so the belt is aligned with the baby's nipples.
- Place the electrodes, Velcro-side down, on either side of the belt as follows:
  - Place the electrode with the **white** lead wire on the baby's **right** side.
  - Place the electrode with the **black** lead wire on the baby's **left** side.
- Place the electrodes far enough apart so that when the belt is wrapped around the baby, the electrode will be located under the arm and lined up with the nipples.
- Be sure the lead wires and patient cable are leading down and away from the baby's face and neck.

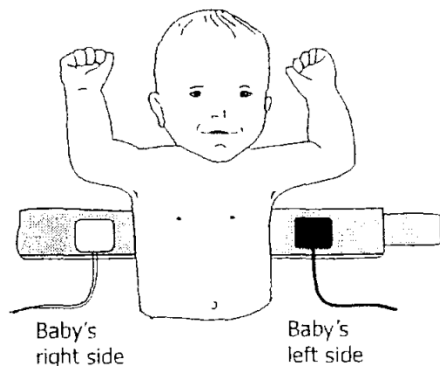


### Step 6: Wrap the electrode belt around the baby.

- Wrap the belt around the baby's chest and fasten it with the Velcro tab.
- The belt should be snug enough so that you can only insert 2 of your fingers (with your hand lying flat against baby) between the belt and the baby.

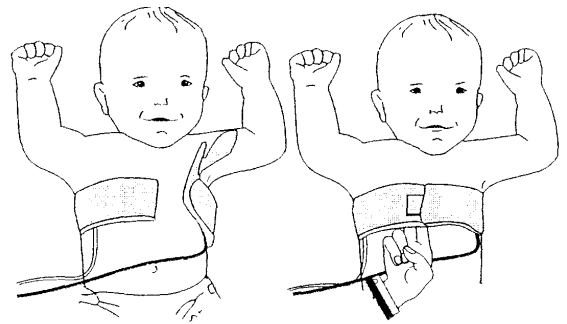

- Insert the **black LA** lead wire into the Velcro side of one electrode. Place the electrode with the **black** lead wire on the baby's **left** side.
- Insert the **white RA** lead wire into the Velcro side of the other electrode. Place the electrode with the **white** lead wire on the baby's **right** side.

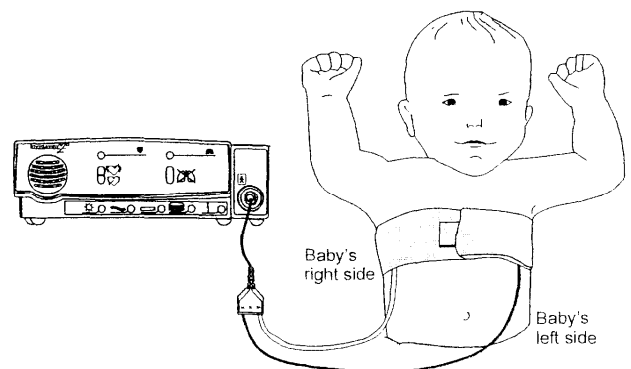


If you have any questions or concerns, please call Cook Children's Home Health at (800) 747-8242. A nurse is available 24 hours a day.

### Disposable Self Adhesive Electrodes

Follow the steps below if you are using disposable electrodes:


- Attach lead wires to the Self Adhesive Electrodes if not already attached.
- Insure the infant's skin is clean and dry.
- Place the electrode with the **white** lead wire on the baby's **right** side, under the arm, 2 finger widths below or lined up with the nipples.
- Place the electrode with the **black** lead wire on the baby's **left** side, under the arm, 2 finger widths below or lined up with the nipples.
- An electrode belt is not needed when using disposable electrodes.

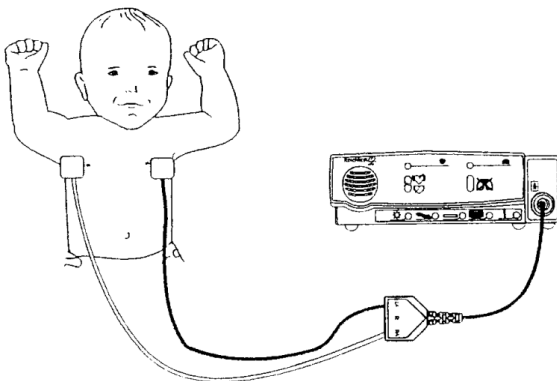

**IMPORTANT:** Do **not** use oils, lotions, or powders on the area of skin that the electrodes will be placed.

### Step 7: Connect the Power Cord/Battery Charger.

- Insert the round connector of the power cord/battery charger into the socket on the back of the panel of the monitor.
- Line up the notch on the connector.
- Push until the connector is fully inserted.
- Plug the power cord/battery charger into a power outlet. The green charge light on the monitor will now come on.

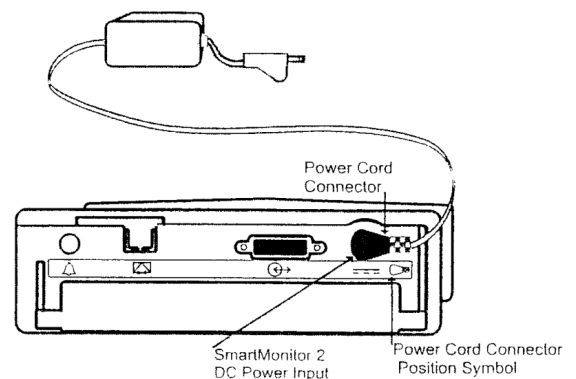

**IMPORTANT:** Do not use monitor if power cord is damaged. Contact Cook Children's Home Health to get the power cord replaced.

If you have any questions or concerns, please call Cook Children's Home Health at (800) 747-8242. A nurse is available 24 hours a day.